

## Benefits of the Painless IPL

- One (1) year warranty
- Canadian Company
- One-on-one training included
- 300,000 pulses per handle
- \$3,500 handle cost is lowest consumable cost compared to competitors (\$15,000 plus)
- Both HR and SR handles are included with the device
- Our IPL handles have a sapphire crystal which distributes heat evenly at every point on the crystal which means overlapping is not necessary. Quicker treatments = more appointment times.
- The largest crystal on the market measuring at 15mm x 50mm
- We have the fastest cooling system out there thanks to our four copper radiators and fans. This enables the Painless to run for eight consecutive hours!
- We have a large, user friendly, multi-language LCD touchscreen

## Can Help With:

- Broken capillaries
- Facial lines and wrinkles
- Collagen stimulation
- Large skin pores
- Hair removal
- Pigmentation
- Skin tightening

MEDSKIN  
DISTRIBUTION



# Technical Specifications of the Painless IPL

- Pulse Mode - single pulse & rapid mode
- Spectrum Range; Skin Rejuvenation 560nm-1200nm  
Hair Removal-690-1200nm

- 1) Technology: Square-Motion Technology
- 2) Handpieces: Two hand pieces  
HR Hand piece: 15\*50mm spot size, 640-1200nm filter  
SR Hand piece: 8\*40mm spot size, 530-1200nm filter
- 3) Display: 10.4 inch colorful touch screen
- 4) Operation Mode: Stationary mode and Motion mode
- 5) Power Supply: 1200W
- 6) Customer Management: Not Available
- 7) Capacitor Capacity: 3 \*10, 000uf
- 8) Repetition Frequency: 2Hz for SR, 3Hz f  
or HR
- 9) Timers: 1S, 3S, 12S for SR  
1S, 3S, 30S for HR.
- 10) Cooling: 1-5 level
- 11) Fluence: 1-7J/cm<sup>2</sup> for HR  
1-12J/cm<sup>2</sup> for SR
- 12) Cooling System: Water + air + semi-conductor
- 13) Dimension: 420\*420\*920mm (with wheels)
- 14) Voltage: 100-120V, 50/60Hz
- 15) Working Temperature: 25°C; Relative humidity: 55%
- 16) Input current: Momentary current 30A/Continuous  
current 4A
- 17) Net Weight: 42kg



604-716-8500  
[www.medskinsupply.com](http://www.medskinsupply.com)


MEDSKIN



**The Painless IPL and how it helps skin rejuvenation** The signs of photo-aging and sun damage affect both men and women and start to appear on the skin from age thirty onward. Most skin changes are commonly associated with aging resulting from sun exposure. These include changes in skin pigmentation such as freckles, melasma, uneven skin tone, rough skin texture, fine lines, wrinkles and spider veins. The signs of sun damage can be treated and reduced with a series of Intense Pulsed Light (IPL) treatment. IPL is very effective in reducing the visible signs of sun damage by stimulating skin renewal, growth of new collagen and elastin fibers in the deeper layers of the skin. This contributes to slowing down the future process of aging. IPL uses the principle of neocollagenesis. Neocollagenesis is a process in which new collagen fibers are created under the skin thus replacing the damaged ones and filling the spaces where the fibers were destroyed, mainly due to environmental factors. IPL treatments enable this process, letting the skin heal by itself while reducing the wrinkles. At the same time, pigmented and vascular lesions are reduced by using a low-energy and broadband light thus the skin imperfections is replaced by newer and younger looking skin. Multiple treatments are required in order to achieve optimum results.

**Collagen Rebuilding/Fine Line Reduction/Scars/  
Dull Skin Rejuvenation**

How it works: IPL therapy destroys undesirable collagen structures. This allows the body to initiate the healing process, which in turn allows it to boost collagen production as well as encouraging blood flow to treated areas.

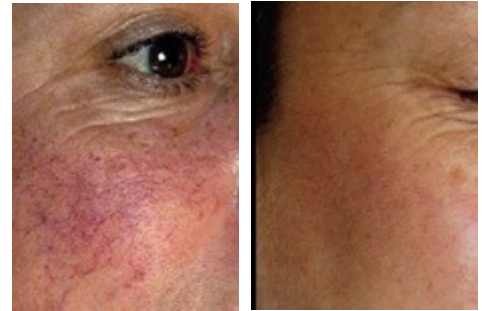

Broken Capillaries courtesy of La Combe Medi-Spa

**Freckle Reduction/Brown Spot Reduction/Age Spot Reduction/  
Moles**

How it works: IPL therapy targets the chromophore through heat and light coagulating the melanin, which eventually sloughs off the skin.




Age Spot Reduction courtesy of Refresh Medi-Spa

**Roseacea/Broken Capillaries/Spider Vein**

How it works: IPL therapy targets the blood vessel destroying it through heat and light. The vessel is dissolved and absorbed naturally into the body.

**Acne**

How it works: IPL therapy targets and kills the bacterial within the pores as well as reducing inflammation associated with lesions.

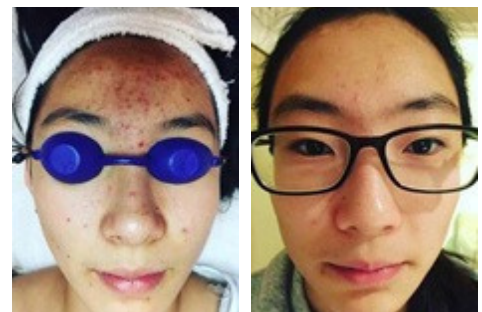

Acne Reduction courtesy of Unwind Spa

If you just did the following Hair Removal and Skin Rejuvenation treatments per month, you could make up to \$10,000 monthly

# of Treatments	Treatment Area	# of Flashes/ Treatment	Total Flashes	Charge/ Treatment	Gross Income	Gross Cost	Net Income
6	Upper Lip	3	18	50	300	3.24	296.76
6	Chin	5	30	65	390	5.40	384.60
2	Full Face	25	50	150	300	9.00	291.00
6	Axilla (Underarms)	10	60	75	450	10.80	439.20
4	Bikini	16	64	90	360	11.52	348.48
4	Brazilian	30	120	120	480	21.60	458.40
4	Upper Legs	200	800	250	1000	144.00	856.00
4	Lower Legs	150	600	250	1000	108.00	892.00
4	Full Legs	350	1400	450	1800	252.00	1548.00
2	Back	200	400	350	700	72.00	628.00
3	Chest	200	600	300	900	108.00	792.00
4	Skin Rejuvenation	200	800	200	800	144.00	656.00
4	Rosacea/Spider Vein	50	200	200	800	36.00	764.00
4	Melasma/Pigmentation	50	200	200	800	36.00	764.00
57			5342		10080	961.56	9118.44

Opportunity  
awaits

For further information call us today at 1-604-716-8500

**MEDSKIN**  
DISTRIBUTION






## Safety, comfort and efficiency with the Painless IPL

The touchscreen interface is pre-programmed for ease of use for the new technician. It is as easy as 1-2-3! As it has 2 hand piece ports, one for Hair Removal and one for Skin Rejuvenation, you can easily switch from one treatment to the other without removing any hand piece. With the Painless IPL, you don't need to be a medical aesthetician to safely and effectively treat the following conditions:


- *Hair Removal*
- *Folliculitis (Ingrown hairs)*
- *Collagen Remodelling*
- *Fine Lines*
- *Dull Skin*
- *Pigmentation*
- *Freckles*
- *Age Spots*
- *Epidermal Melasma*
- *Rosacea*
- *Broken Capillaries*
- *Scars*

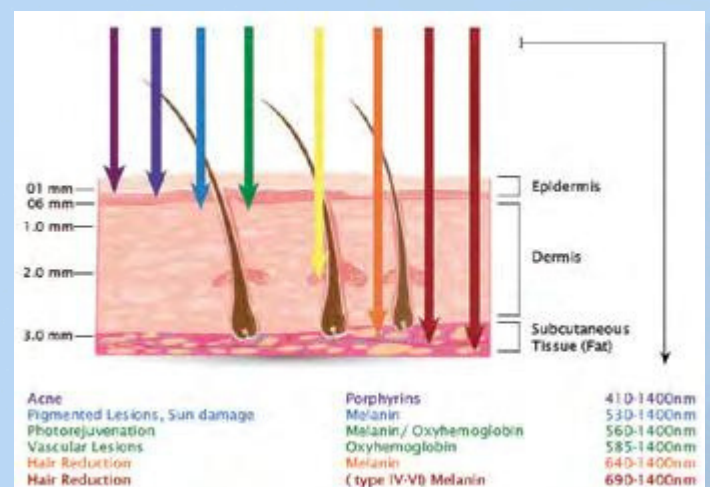

The Painless IPL large sapphire crystal head allows you to treat your clients quickly and the built in cooling system keeps them comfortable.

The size of the Painless Hair Removal hand piece is 15mm x 50 mm making it one of the largest in the market. The spot size of the Skin Rejuvenation hand piece is 8mm x 40mm, a perfect size to treat all of the intricacies on the face. At just \$3,500 per handle and up to a 300,000 pulse life, it has the lowest consumable cost compared to competitors.

## HOW LASER HAIR REMOVAL WORKS

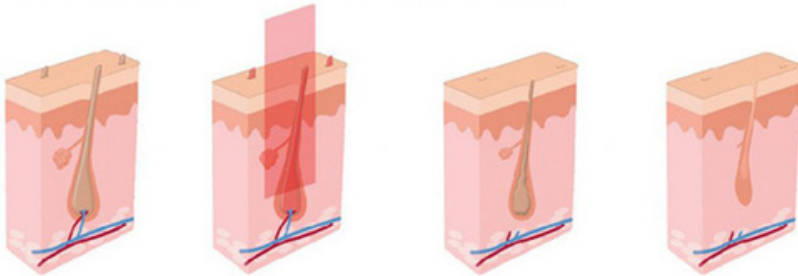

1 at your initial treatment the hair should be visible on the surface of the skin.

2 As the laser is directed onto the hair heat energy is drawn into the melanin, which in turn heats the rest of the hair shaft.

3 The heat either weakens or completely destroys the hair follicle. As a result the follicle detaches from the hair shaft.

4 The hair sheds and, without the follicle, hair cannot regrow in the area.

Actual size of the Painless IPL  
Hair Removal crystal

Fast, effective  
and effortless

